



PERSON LIFT HANDYLIFT

COD. SPGH-MDE

Handylift is an electro-mechanical lifting arm which facilitates the passengers transfer from their wheelchair to the original vehicle seat. It is normally fitted in line with the front passenger door of the vehicle, however it can be positioned in other entry points of vehicles which are suitable for its installation.

The harness, made of breathable and washable material, is hooked to the **Handylift** which, with its three horizontal rotation points, and electronic lifting system, allows to the person to be transferred quickly and easily.

The latest version of this device has a weight of 9,5 kg, a lifting capacity of 120 kg (reduced to 90 kg if used with hinge door bracket), and has a WiFi remote control to operate the electric motor. It is equipped with an anti-crush system and LED lighting to facilitate loading and unloading.

Handylift is manufactured with treated and powder coated steel, and the rotating bushings are made from a techopolymer which, along with the IP66 rated motor and IP44 electronic unit, allows it to be used in damp or humid environments.

Once the transfer is complete, **Handylift** can be removed from the support bracket fitted to the vehicle and placed in its own bag. **Handylift**, being a single device, allows transfers to be carried out on different vehicles which have been fitted with a hinge door bracket.

The **Handylift** device, now officially recognized as Class I Medical Device (Conform to 93/42 CE directive), can be used in combination with different optional, specific made for home use.

PERSON LIFT HANDYLIFT

COD. SPGH-MDE



- Equipped with fixed or removable bracket.
- All detachable.
- Supplied with its own comfortable box.
- Washable and ergonomic harness available in three sizes (S - M - L).
- Capacity → 120 Kg (90 KG with hinge door bracket).
- Lifting excursion → 440 mm.
- Electric movement with Wi-Fi remote control and manual emergency system.
- LED lighting.



WHY SPGH-MDE ?

- ✓ **No crushing system**
- ✓ **Wi-Fi remote control**
- ✓ **Installable on every vehicle**
- ✓ **LED lighting for nights use**
- ✓ **Weight only 9,5 KG**





Marchio di **CARROZZERIA 71**srl
Strada Patro, 3 - 14036 Moncalvo (AT) Italia
Tel. +39 0141 91 76 65 - Fax +39 0141 92 39 79
info@handytech.it - www.handytech.it




Member of

Approved by

Production process certified by

Certified company

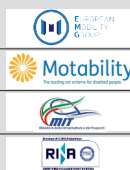

Dealer

